


Performance Verification



Eliminate error-prone manual data transcription and conduct statistical analysis to deliver performance verification reports

ACL TOP Family 70 Series Hemostasis Testing Systems, together with HemoHub Intelligent Data Manager, deliver automated performance verification tools that simplify study development and execution workflow, eliminating manual steps.

Transform workflow from multiple manual steps

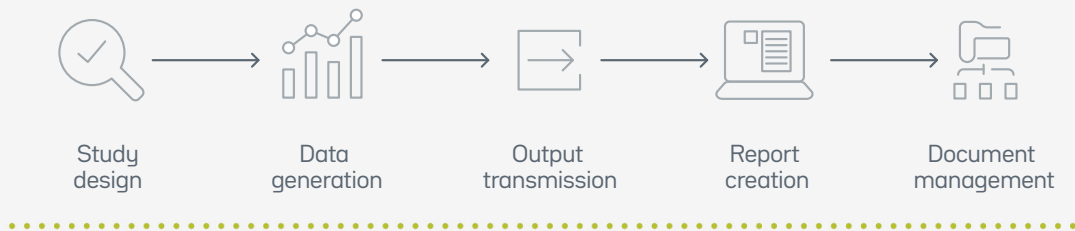

To just two automated steps

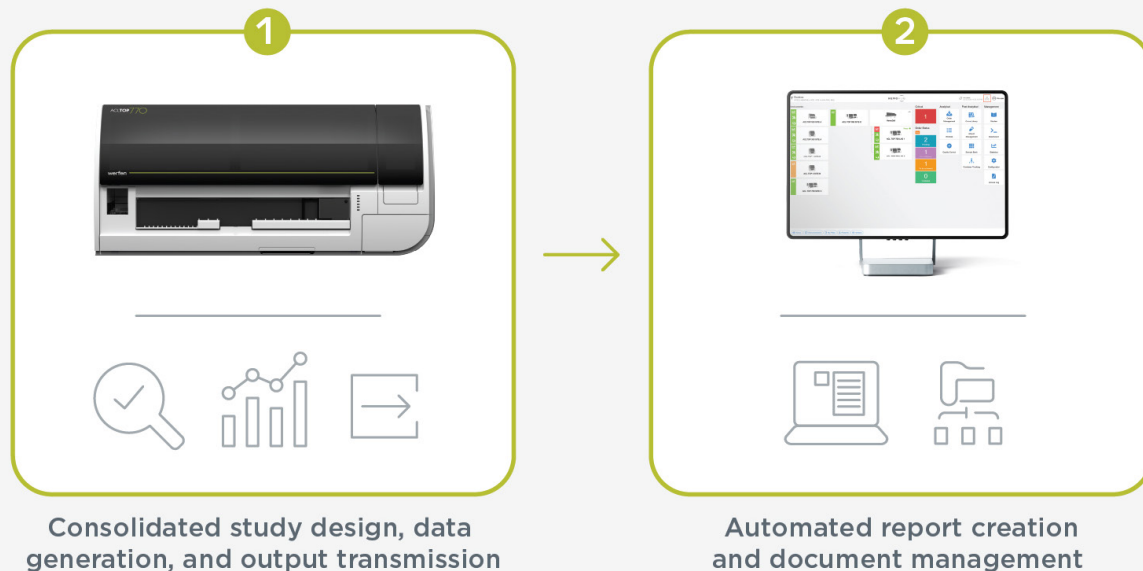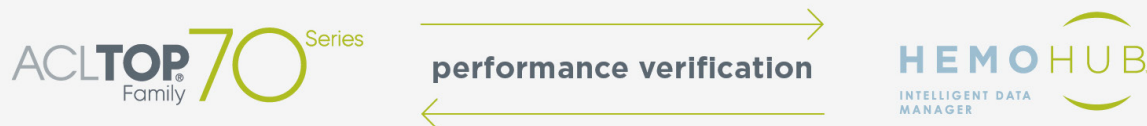

Automate key studies

Seamless automated workflows for key studies include:

Reagent Lot Conversion

- Perform reagent lot conversion studies with multiple lots onboard simultaneously
- In safe data collection mode, run samples for both active and alternate lots, without risk of transmitting study data to lab information system (LIS)

Reference Interval

- Calculate Mean of Normal (reference range) with option to apply across participating ACL TOP Family 70 Series instruments (by authorized personnel), with the click of a button


Instrument Comparison

- Compare networked ACL TOP Family 70 Series instruments, connected to HemoHub. Following study creation, results are compared within HemoHub, and the corresponding report can be transmitted to every participating instrument



Automate report creation and document management with HemoHub

- Enables acceptance-criteria setting, design and layout of report content, and automation of final study reports
- Final reports with study conclusions are automatically stored on HemoHub and instrument for quick access

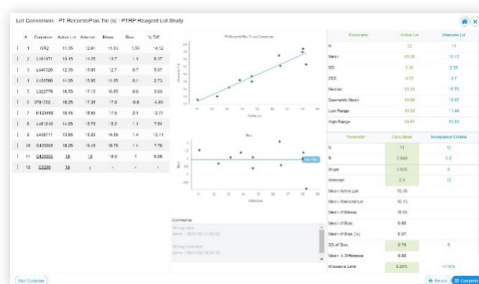

Fast and flexible for enhanced lab efficiency.

werfen.com

For more information, contact your local Werfen sales representative or distributor.

ACL, ACL AcuStar, ACL Elite, ACL TOP, HemoCell, HemoHub, HemosIL, ReadiPlasTin, RecombiPlasTin, SynthAFax, and SynthASil are trademarks of Instrumentation Laboratory Company (d.b.a. Werfen) and/or one of its subsidiaries or parent companies and may be registered in the United States Patent and Trademark Office and in other jurisdictions. The Werfen logo is a trademark of Werfen and may be registered in the Patent and Trademark Offices of jurisdictions throughout the world. All other product names, company names, marks, logos, and symbols are trademarks of their respective owners.

©2026 Instrumentation Laboratory. All rights reserved.

werfen